

Pd(II) EnCat™

Encapsulated Palladium(II) Catalysts

Reaxa's Pd EnCat™ catalysts immobilise palladium acetate together with activating ligands within a robust porous polyurea matrix to allow easier, faster and cleaner processes to be developed

Contaminated crude Suzuki coupling product made with homogeneous Pd acetate catalyst system



Clean product made using Reaxa's Pd EnCat™ 30 catalyst with no extra purification

Cleaner products

Cleaner waste streams

Fast, efficient processes

No plant contamination

Improved processes

Process intensification

typically less than 10 ppm Pd in crude product

reduced metal losses in Pd EnCat™ processes

the EnCat™ beads filter easily

metal & ligands are trapped within the beads


high activity in many types of coupling reactions

can be used in batch & flow processes


<i>Product</i>	<i>Pd content % w/w</i>	<i>Co-encapsulated ligand</i>
Pd(II) EnCat™ 30	4.3	none
Pd(II) EnCat™ 40	4.6	none
Pd(II) EnCat™ TPP30	4.7	
Pd(II) EnCat™ TOTP30	4.7	
Pd(II) EnCat™ BINAP30	4.7	
Pd(II) EnCat™ polyTPP30	4.6	

Pd EnCat™ Applications

Pd EnCat™ Chemistry Examples:

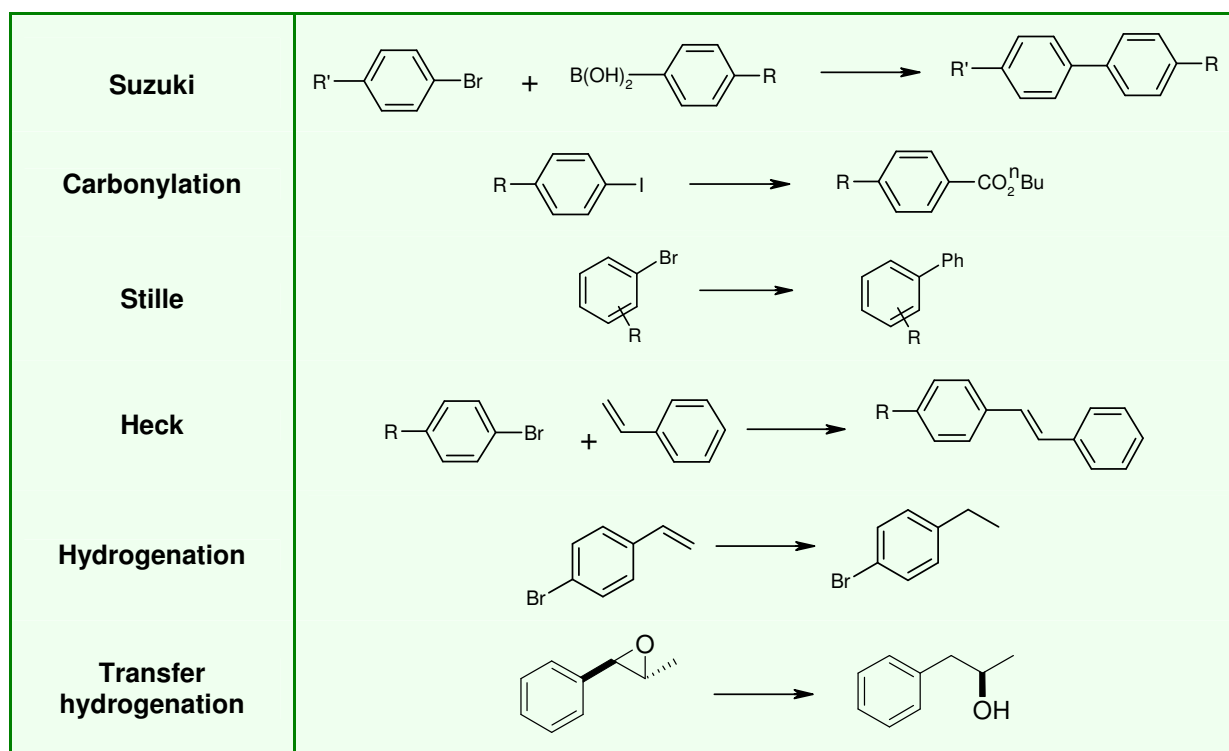

Process Scale Suzuki Coupling Example:

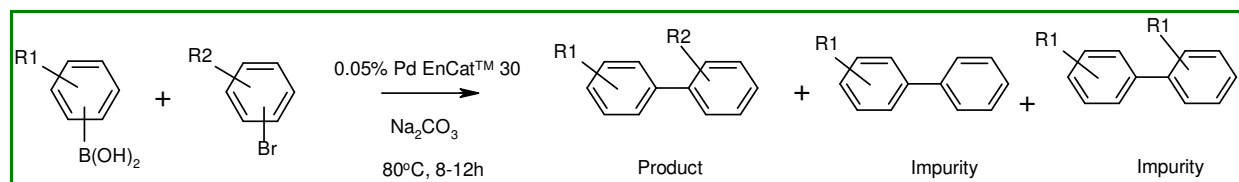

The client's process using 2.5 mol% Pd/C gave a 94% yield with impurities at 5% and 1% and a reaction time of 60 mins. Following extensive work-up and recrystallisation, Pd levels of > 50 ppm were present in the product. A switch to Pd(II) EnCat™ 30 at 0.25 mol% gave a dramatic improvement: after 30 mins the crude product was reproducibly obtained at >99% yield, with both impurities reduced to <1% and with Pd levels reduced to <10 ppm without recrystallisation, representing a significant cost saving for the client.

Selected References:

1. C.K.Y. Lee, Connie, K.Y. Lee, A.B. Holmes, S.V. Ley, I.F. McConvey, B. Al-Duri, G.A. Leeke, R.C.D. Santos, J.P.K. Seville, *Chemical Communications*, 2005, (16), 2175 – 2177.
2. I.R. Baxendale, C.M. Griffiths-Jones, S.V. Ley and G. Tranmer, *Chem. Eur. J.*, 2006, 12, 4407-4416.
3. J. Gil-Molto, S. Karlstrom and C. Najera, *Tetrahedron*, 2005, 61, 12168–76.

R&D quantities of Pd EnCat™ products are available from Sigma-Aldrich at:
www.sigma-aldrich.com

For technical support, bulk quotations & information on Pd EnCat™ catalysts please contact:
info@rexa.com

Rexa Limited, Hexagon Tower, Blackley, Manchester, M9 8ZS, UK

Tel: +44 (0)161 721 2202 Fax: +44 (0)161 721 2005 www.rexa.com

Company Registered in England No. 04884654